Orthofix US LLC  
 3451 Plano Parkway  
 Lewisville, Texas 75056-9453 U.S.A.  
 1-214-937-3199  
 1-888-298-5700  
 www.orthofix.com

Australian Sponsor  
 Emergo Australia  
 Level 20, Tower II  
 Darling Park  
 201 Sussex Street  
 Sydney, NSW 2000  
 Australia

Medical Device Safety  
 Service (MDSS)  
 Schiffgraben 41  
 30175 Hannover  
 Germany  
 www.mdss.com

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

Click directory below for desired language

|            |    |         |
|------------|----|---------|
| English    | en | 2-13    |
| Español    | es | 14-27   |
| Português  | pt | 28-41   |
| Français   | fr | 42-55   |
| Italiano   | it | 56-68   |
| Deutsch    | de | 69-82   |
| Česky      | cz | 83-95   |
| Nederlands | nl | 96-108  |
| Magyar     | hu | 109-121 |
| Ελληνικά   | el | 122-136 |

**INSTRUCTIONS FOR USE**  
**Important Information – Please Read Prior to Use**

**Rx Only**  
**CE** 2797

Orthofix US LLC  
3451 Plano Parkway  
Lewisville, Texas 75056-9453 U.S.A.  
1-214-937-3199  
1-888-298-5700  
www.orthofix.com

Medical Device Safety Service  
(MDSS)  
Schiffgraben 41  
30175 Hannover  
Germany  
+49 511 6262 8630  
www.mdss.com

**English** **EN**

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

### Description:

The Orthofix surgical instruments covered by this insert are designed for use in orthopedic and spine surgical procedures for implantation of devices manufactured by Orthofix US LLC. Prior to use of these instruments, please refer to the product instructions for use specific to the product line and the surgical procedure to be followed.

**Note:** Personnel responsible for the reprocessing of reusable Instruments shall have appropriate training and competency in accordance with the Health Care Facility's procedures or seek out additional education as needed.

### Warnings and Precautions:

- Do not reprocess single-use surgical instruments.
- Implants that are provided non-sterile may be sterilized multiple times before their first clinical use
- Implants are not designed to be cleaned or disinfected; any implant that comes in contact with bodily tissue or fluid or otherwise becomes soiled must be disposed of in accordance with applicable policies and regulations
- Reprocess the reusable surgical instruments as soon as is reasonably practicable in order to minimize the drying of soil and residuals. For optimal results, instruments should be cleaned within 30 minutes of use.
- Use a washer-disinfector in compliance with EN ISO 15883 that is properly installed, qualified and regularly subjected to maintenance and testing.
- Ensure that the cleaning receptacle is clean and dry, no visible foreign material shall be present prior to use.
- Ensure that the washer-disinfector and all services are operational prior to use.
- Personnel that work with contaminated medical devices must follow safety precautions as per the procedure of the healthcare facility.
- Cleaning and disinfection solutions with a pH 7-10.5 are recommended.
- Cleaning and disinfection solutions with a higher pH should be used in accordance to the material compatibility requirements stated on the detergent technical data sheet.
- Aluminum based instruments are damaged by alkaline (pH > 7) detergents and solutions.
- Detergents and disinfectants with fluoride, chloride, bromide, iodide or hydroxyl ions must not be used.
- Contact with saline solutions should be minimized.
- Certain cleaning solutions such as those containing caustic soda, formalin, glutaraldehyde, bleach, and/or other alkaline cleaners may damage instruments. These solutions should not be used.
- Complex devices such as those with hinges, lumens or mated surfaces must be thoroughly manually pre-cleaned before automated washing in order to remove soiling that accumulates in recesses.
- Do not use metal brushes or steel wool.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

- Carefully inspect all reprocessed instruments prior to use for possible damage, unacceptable wear, or non-functioning components.
- Instruments must be cleaned and disinfected before they are returned to the manufacturer for any reason.

### **Introduction:**

This document gives general guidance on how medical devices supplied by Orthofix may be reprocessed to prepare them for use.

The efficacy of the reprocessing also depends on the equipment, operators, cleaning agents and procedures used.

It is the responsibility of the health care facility to ensure that the reprocessing is performed in accordance with the instructions provided. Alternative methods of reprocessing are allowed, but it is the responsibility of the health care facility to validate and monitor the process. Any deviation from the instructions provided by Orthofix must be validated by the end user.

Prior to use and reprocessing, Orthofix recommends that users take into account any international, national or local requirements regarding reusable devices and their reprocessing. This includes any requirements related to equipment (e.g. steam autoclaves, washer disinfectors, detergents, disinfectants, other accessories) and water quality.

The reprocessing instructions reported in the present document have been validated by Orthofix and they are in accordance with the following standards and Technical Information Report:

- |               |                |               |
|---------------|----------------|---------------|
| • AAMI TIR12  | • EN 868-2     | • ISO 17664-1 |
| • AAMI TIR55  | • ISO 11138-7  | • ISO 17665-1 |
| • AAMI TIR68  | • ISO 10993-5  | • ISO 15883-1 |
| • AAMI ST79   | • ISO 10993-17 | • ISO 15883-2 |
| • AAMI ST-98  | • ISO 10993-18 | • ISO 15883-5 |
| • AAMI ST-108 | • ISO 11607-2  | • ISO 21726   |


### **Limitations on Reprocessing:**

- *Repeated reprocessing has minimal effect on reusable instruments.*
- *End of life is commonly determined by worn and/or damaged components, due to use.*

### **Reprocessing Instructions Overview:**

Steps to prepare a reusable medical device for re-use are summarized in the following chart. More detailed instructions for each cleaning step are given below.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide
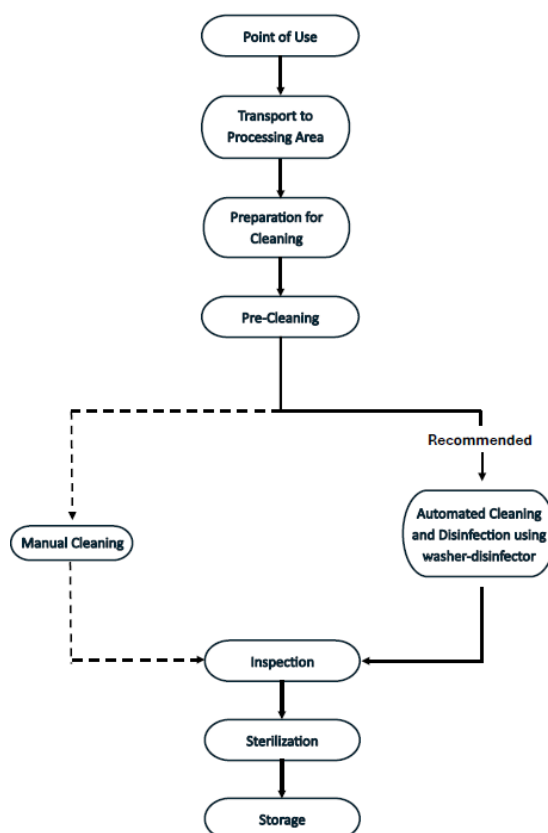
### **Cleaning and Disinfection:**

All reuseable surgical instruments must be disinfected and cleaned using established methods before sterilization and reintroduction into the sterile surgical field. The following recommendations are for manual cleaning and automated cleaning and disinfection of surgical instruments. These recommendations are based on validated laboratory testing; any deviations must be validated by user.

The health care facility is responsible for ensuring that conditions essential to safe handling and decontamination can be achieved.

ANSI/ AAMI ST35 Safe Handling and Biological Decontamination of Reusable Medical Devices in Health Care Facilities and in Nonclinical Settings and ISO 15883-1 Washer-disinfectors – Part 1: General requirements, terms and definitions and tests, provide relevant guidelines.

### **General Considerations:**

Prior to Cleaning, refer to the system Instructions for Use for disassembly instructions for any instrument capable of being disassembled.

Orthofix provides two methods of cleaning: a manual method and an automated method.

Wherever applicable, the cleaning phase should start immediately after the pre-cleaning phase to avoid soil drying.

Staff shall use protective equipment following the safety precautions to comply with the procedure of the healthcare facility. In particular, staff should take note of the instructions provided by the cleaning agent manufacturer for correct handling and use of the product.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

Observe all instructions provided by the detergent's manufacturer regarding the immersion time of the device in the cleaning agent and its concentration.

The quality of the water used for diluting cleaning agents and for rinsing medical devices shall be carefully considered.

### Non-European Cleaning Procedures

#### From Point of Use:

Whenever possible, do not allow blood, debris or body fluids to dry on instruments. For best results and to prolong the life of the surgical instrument reprocess immediately after use.

1. Remove excess body fluids and tissue from instruments with a disposable, non-shedding wipe. Place instruments in a basin of purified water or in a tray covered with damp towels. Do not allow saline, blood, body fluids, tissue, bone fragments or other organic debris to dry on instruments prior to cleaning.
2. For optimal results, instruments should be cleaned within 30 minutes of use or after removal from solution to minimize the potential for drying prior to cleaning.
3. Used instruments must be transported to the central supply in closed or covered containers to prevent unnecessary contamination risk.

**Note:** Soaking in proteolytic enzymatic detergents or other pre-cleaning solutions facilitates cleaning, especially in instruments with complex features and hard-to-reach areas (e.g. cannulated and tubular designs, etc.). These enzymatic detergents as well as enzymatic foam sprays break down protein matter and prevent blood and protein based materials from drying on instruments. Manufacturer's instructions for preparation and use of these solutions should be explicitly followed.

#### Preparation for Cleaning:

1. All instruments with moving parts (e.g., knobs, triggers, hinges) should be placed in the open position to allow access of the cleaning fluid to areas that are difficult to clean.
2. Soak the instruments for a minimum of 10 minutes in purified water prior to the manual or automated cleaning process.
3. Use a soft cloth or a soft plastic bristle brush to remove any visible soil from the instruments prior to manual or automated cleaning. Use a soft plastic bristle brush or a pipe cleaner to remove soil from any inner lumens. You can also use a syringe (if appropriate) for hard to reach areas.
4. Enzymatic detergent should be used for manual and automated cleaning. All enzymatic detergents should be prepared at the use dilution and temperature recommended by the manufacturer. Softened utility water may be used to prepare the enzymatic detergents. Use of recommended temperatures is important for optimal performance of enzymatic detergent.

#### Manual Cleaning:

1. Completely submerge instruments in an enzymatic detergent and allow to soak for 20 minutes. Use a soft-bristled, nylon brush to gently scrub the device until all visible soil has been removed. Particular attention must be given to crevices, lumens, mated surfaces, connectors and other hard-to-clean areas. Lumens should be cleaned with a long, narrow, soft-bristled brush (i.e. pipe cleaner brush).
2. Remove the instruments from the enzymatic detergent and rinse in utility water for a minimum of 3 minutes. Thoroughly and aggressively, flush lumens, holes and other difficult to reach areas.
3. Place prepared cleaning solution in a sonication unit. Completely submerge device in cleaning solution and sonicate for 10 minutes.
4. Rinse instrument in purified water for at least 3 minutes or until there is no sign of blood or soil on the device or in the rinse stream. Thoroughly and aggressively, flush lumens, holes and other difficult to reach areas.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

5. Repeat the sonication and rinse steps above.
6. Remove excess moisture from the instrument with a clean, absorbent and non-shedding wipe.
7. Inspect the instruments for visible soil.
8. If visible soil is noted, repeat the steps listed above.

### Automated cleaning:

1. Completely submerge the instruments in an enzymatic detergent and allow to soak and sonicate for 10 minutes each. Use a soft nylon bristled brush to gently scrub the device until all visible soil has been removed. Particular attention must be given to crevices, lumens, mated surfaces, connectors and other hard to clean areas. Lumens should be cleaned with a long, narrow, soft nylon bristled brush (i.e. pipe cleaner). Use of a syringe or water jet will improve flushing of difficult to reach areas and closely mated surfaces.
2. Remove instruments from the cleaning solution and rinse in purified water for a minimum of 1 minute. Thoroughly and aggressively flush lumens, blind holes and other difficult to reach areas.
3. Place instruments in a suitable washer/disinfector basket and process through a standard instrument washer/disinfector cleaning cycle.
4. Orient instruments into the automated washer's carriers as recommended by the washer manufacturer.
5. The following minimum parameters are essential for thorough cleaning.
  - 5.1. 2 minute prewash with cold utility water
  - 5.2. 1 minute prewash with hot utility water
  - 5.3. 2 minutes detergent wash with hot utility water (64-66°C/146-150°F)
  - 5.4. 1 minute hot utility water rinse
  - 5.5. 2 minute thermal rinse with purified water (80-93°C/176-200°F)
  - 5.6. 1 minute purified water rinse (64-66°C/146-150°F)
  - 5.7. 7 to 30 minute hot air dry (116°C/240°F)
6. Inspect the instruments for visible soil.
7. If visible soil is noted, repeat the above listed steps until no visible soil is noted.

## European Cleaning and Disinfection Procedures

### From Point of Use:

Whenever possible, do not allow blood, debris, or body fluids to dry on instruments. For best results and to prolong the life of the surgical instrument, reprocess immediately after use.

1. All used surgical instruments must be regarded as contaminated
2. Remove excess body fluids and tissue from instruments with a disposable, non-shedding wipe. Place instruments in a basin of utility water or in a tray covered with damp towels. Do not allow saline, blood, body fluids, tissue, bone fragments or other organic debris to dry on instruments prior to cleaning.
3. For optimal results, instruments should be cleaned within 30 minutes of use or after removal from solution to minimize the potential for drying prior to cleaning.
4. Put the instrument in utility water (< 40°C or < 104°F) or a neutral pH cleaning solution for at least 10 minutes, removing gross soil with soft cloth or soft nonmetallic brush.

**Note:** Do Not use a fixating detergent or water > 104°F because this can cause the fixation of residue.

5. Used instruments must be transported to the central supply in closed or covered containers to prevent unnecessary contamination risk.
6. Follow all hospital protocols for handling and transporting contaminated and bio-hazardous materials.
7. Handling, collection and transportation of used instruments must be strictly controlled to minimize any possible risks to patients, personnel and any area of the healthcare facility.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

### Manual Pre-cleaning:

1. Perform a pre-cleaning due to the contamination of the device. Take special care when the items have:
  - Cannulations
  - Long blind holes
  - Mating surfaces
  - Threaded components
  - Rough surfaces
2. Wear protective equipment and follow the safety precautions of the healthcare facility.
3. Ensure that the cleaning receptacle is clean and dry, and no visible foreign material can be present.
4. Fill the receptacle with sufficient detergent solution.  
Orthofix recommends the use of a slightly alkaline enzymatic detergent solution based on a detergent containing <5% anionic surfactants and enzymes, prepared using deionized water, at room temperature (18° C-25° C) (e.g. 2% concentration, Neodisher MediZyme).
5. Carefully immerse the component in the solution in order to displace trapped air.
6. Scrub the device in the cleaning solution with a soft bristled nylon brush until all visible soiling is removed. Use a soft bristled nylon brush to remove residuals from lumens, rough or complex surfaces using a twisting motion.
7. Rinse cannulations with cleaning solution using a syringe. Never use metal brushes or steel wool.
8. Remove the device from the cleaning solution.
9. Brush the single components in utility water.
10. Ultrasonic clean the single components in a degassed cleaning solution. Orthofix recommends, on the basis of the validation performed, the use an ultrasound frequency of 25-40 kHz for at least 4 minutes.
11. Rinse the components in critical water until all traces of cleaning solution are removed. Use a syringe in case of lumens or cannulations.
12. Remove item from rinse water and drain.
13. Carefully hand-dry using clean, lint free cloth.

### Manual Cleaning:

1. Wear protective equipment following safety precautions to comply with the procedure of the healthcare facility.
2. Ensure that the cleaning receptacle is clean and dry, and no visible foreign material can be present.
3. Fill the receptacle with sufficient room temperature (18°C-25° C) cleaning solution.  
Orthofix recommends the use of a slightly alkaline enzymatic cleaning solution (e.g. 2% concentration, Neodisher MediClean forte).
4. Carefully immerse the component in the solution in order to displace trapped air; it is important to ensure that the cleaning solution reaches all surfaces, including holes or cannulations. Conduct the immersion phase for 6-8 minutes.
5. Thoroughly scrub the device in the cleaning solution with a soft bristled nylon brush until all visible soiling is removed. Use a soft bristled nylon brush and a twisting motion to remove residuals from lumens, rough or complex surfaces.
6. Rinse cannulations at least three times with cleaning solution using a syringe. Never use metal brushes or steel wool.
7. Remove the device from the cleaning solution.
8. Brush the single components in utility water.
9. Put single components in ultrasonic device with degassed cleaning solution:
  - 9.1. Submerge in appropriate solution at 40-60°C.  
Orthofix validated the use of a detergent solution based on a detergent containing <5% anionic surfactants, non-ionic surfactants and enzymes, prepared using deionized water (e.g. 2% concentration, Neodisher MediClean forte).

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

- 9.2. Orthofix recommends, on the basis of the validation performed, to use an ultrasound frequency of 25-40 kHz for 15 minutes.

The use of other solutions and parameters shall be validated by the user and the concentration shall be in compliance with the detergent's manufacturer technical datasheet.

10. Rinse the components in critical water until all visible traces of cleaning solution are removed.
11. Rinse the cannulations, rough or complex surfaces at least three times with critical water. When cannulations are present use a syringe to facilitate this step.
12. Remove item from rinse water and drain.
13. If, after completion of the cleaning steps, some encrusted soil remains on the device, the cleaning step must be repeated as described above.
14. Carefully hand-dry using clean, lint-free cloth.

### Automated Cleaning and Disinfection:

1. Use a washer-disinfector in compliance with EN ISO 15883 that is properly installed, qualified and regularly subjected to maintenance and testing
2. Ensure that the cleaning receptacle is clean and dry, and no visible foreign material can be present.
3. Ensure that the washer-disinfector and all services are operational.
4. Load the medical devices into the washer-disinfector.
  - 4.1. Ensure products are in the open position or disassembled, as appropriate, before placing them in the baskets.
  - 4.2. Wherever possible, all parts of disassembled devices should be kept together in one container.
  - 4.3. Place heavier devices in the bottom of the baskets.
5. Connect cannulations to the rinsing ports of the washer-disinfector.
  - 5.1. If no direct connection is possible, locate the cannulations directly on injector jets or in injector sleeves of the injector basket.
  - 5.2. Orient instruments into the automated washer's carriers as recommended by the washer manufacturer.
6. Avoid contact between devices because movement during washing could damage the devices and the washing action could be compromised.
7. Arrange medical devices so that cannulations are in vertical position and blind holes incline downwards to promote the drainage of any material.
8. Select and start a cycle according to the recommendations of the washer manufacturer.
9. Use a wash and thermal disinfection program validated by the washer-disinfector manufacturer. When using alkaline solutions, a neutralizer must be added. Orthofix recommends the following steps at a minimum:
  - 9.1. Pre-cleaning for 4 min.
  - 9.2. Cleaning with the appropriate solution for 10 min at 40-60°C.

Orthofix validated the use of a detergent solution based on a detergent containing <5% anionic surfactants, non-ionic surfactants and enzymes, prepared using deionized water (e.g. 0.5% concentration, Neodisher MediClean forte).
  - 9.3. Neutralize with a neutralizing agent solution.

Orthofix validated the use of a 0.1% citric acid based agent solution for 6 minutes (e.g. 0.1% concentration, Neodisher Z).
  - 9.4. Final rinsing with critical water for 3 minutes.
  - 9.5. Thermal disinfection at minimum 90°C or 194°F (max 95 °C or 203°F) for at least 1 minute or until  $A_0 = 600$  is reached.
  - 9.6. The water used for thermal disinfection must be critical water.
  - 9.7. Drying at 110°C for 40 minutes.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

The suitability of other solutions, concentration, time and temperature shall be checked and validated by the user following the detergent manufacturer's technical datasheet.

**Note:** When the instrument has a cannulation, an injector should be used to dry the internal part.

10. On completion of the cycle, ensure that all stages and parameters have been achieved.
11. Wearing protective equipment, unload the washer disinfector when it completes the cycle.
12. If necessary, drain excessive water and dry by using clean, lint-free cloth.
13. Visually inspect each device for remaining soil and dryness. If soil remains, repeat the cleaning process as described above.

### Instructions for Assembly:

Refer to the system operative technique for assembly instructions for any instrument capable of being disassembled for reprocessing. Disassembly of instruments after functional testing may be required based on sterilization tray configuration.

### Cleaning Agent Information:

Orthofix used the following cleaning agents during validation of these processing recommendations. These cleaning agents are not listed in preference to other available cleaning agents which may perform satisfactorily.

- For manual precleaning: Neodisher Medizym concentration 2%
- For manual cleaning: Neodisher Mediclean forte concentration 2%
- For automated cleaning: Neodisher Mediclean forte concentration 0.5%; Neodisher Z concentration 0.1% for neutralization

### Maintenance, Inspection and Functional Testing:

The following guidelines shall be applied to all Orthofix instruments that are labeled for multiple use. All functional checks and inspections described below also cover the interfaces with other instruments or components to identify and/or prevent failure modes caused by end of life of the product, improper use or improper maintenance.

Orthofix does not typically specify the maximum number of uses for re-useable medical devices. The useful life of these devices depends on many factors, including the method and duration of each use and the handling between uses. Careful inspection and functional testing of the device before use is the best way to determine whether a device has reached the end of its serviceable life.

DO NOT use instruments which have reached the end of their serviceable life. Discard the instrument in accordance with facility procedures and any applicable policies and regulations or return to Orthofix for disposal.

As a general preventive action Orthofix recommends consulting the operative technique for any special use or care instructions to avoid damage related to incorrect use.

The following general instructions apply to all Orthofix products:

- All instruments and product components must be visually inspected under good light for cleanliness.
- If visual inspection indicates that the device was not properly cleaned, repeat the cleaning and disinfection steps.
- If the device is still not properly cleaned after repeated attempts, by visual inspection, discard the device in accordance with facility procedures and any applicable policies and regulations or return to Orthofix for disposal.
- All instruments and product components must be visually inspected for any signs of deterioration that may cause failure during use (such as rust, pitting, scratches, or cracks), and its functions must be tested

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

before being sterilized. If a component or instrument is believed to be faulty, damaged or suspect, it must NOT BE USED.

- Products that show excessive fading of marked product code, UDI and lot, thus preventing clear identification and traceability, must NOT BE USED.
- Instruments intended for cutting bone and/or tissue (e.g. tap, rasp, curette, rongeur) must be checked for sharpness and when any of the cutting surfaces show signs of wear such as nicks, abrasions or otherwise dulled cutting surfaces must NOT BE USED.
- All calibrated instruments must be monitored and maintained according to their respective calibration schedules
- When instruments form part of an assembly, check assembly of components and designated locations within their respective sterilization tray
- All instruments must be placed back in their proper tray locations before packaging and sterilization.

It is important to follow the cleaning procedure suggested by Orthofix to avoid damage related to incorrect handling.

### Torque Limiting Instrument Maintenance

- If a torque-limiting handle has been dropped, impacted or used incorrectly, return to Orthofix.
- Torque-limiting handles require maintenance at minimum, every three years or per your service agreement. Please return your torque limiting handles to Orthofix for required maintenance.

### **Sterilization:**

#### **Disclaimer:**

Orthofix instrument cases and trays are intended to protect instrumentation and facilitate the sterilization process by allowing steam penetration and drying. Orthofix has verified through laboratory testing that its instrument cases and trays are suitable for the specific sterilization methods and cycles for which they have been tested. Health care personnel bear the ultimate responsibility for ensuring that any packaging method or material, including a reusable rigid container system and sterilization wraps, is suitable for use in sterilization processing and sterility maintenance in a particular health care facility. Testing should be conducted in the health care facility to ensure that conditions essential to sterilization can be achieved. Orthofix does not accept responsibility or liability arising from a lack of cleanliness or sterility of any medical devices supplied by Orthofix that should have been properly cleaned and/or sterilized by the end user prior to use.

#### **Responsibility of the user:**

Non-sterile surgical instruments and implants are to be sterilized by the hospital prior to surgery. These instructions are not applicable to instruments or implants that are provided sterile.

Health care personnel bear the ultimate responsibility for ensuring that any packaging method or material is suitable for use in sterilization processing and sterility maintenance. Please consult local regulations and comply with any specific rules concerning the reprocessing of reusable devices.

Users should conduct testing in the health care facility to ensure that conditions essential to sterilization can be achieved and that specific configuration of the container contents is acceptable for the sterilization process and for the requirements at the point of use.

ANSI/AAMI ST33, *Guidelines for the Selection and Use of Reusable Rigid Container Systems for Ethylene Oxide Sterilization and Steam Sterilization in Health Care Facilities*, and ISO 17665-1, *Sterilization of health care products – Moist heat – Part 1: Requirements for the development, validation and routine control of a sterilization process for medical devices*, provide relevant guidelines.

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

### Packaging and Sterilization Process:

- Implants and instruments should be placed in the designated location in their respective tray.
- Trays should be double wrapped then sterilized to maintain sterility (see below for recommended wrapping configuration).
- Gravity cycle is not recommended in Europe, it is not considered state of the art.
- Cycles at 132°C are not usually available in Europe, therefore please use only available validated cycles.

Trays are to be double wrapped using FDA-cleared sterilization wraps (2 wraps) prior to placement in autoclave. In Europe, an EN 868-2 compliant sterilization wrap may be used; folded to create a sterile barrier in accordance to a process validated per ISO 11607-2. The recommended sterilization cycle will produce a Sterility Assurance Level (SAL) of at least  $10^{-6}$ .

Orthofix recommends always using a pre-vacuum cycle for steam sterilization. The Gravity cycle was validated but is suggested only when no other options are available. The gravity cycle was only validated for wraps and was not validated for sterilization in rigid containers.

In a properly functioning and maintained steam sterilizer (per ISO 17665-1 guidelines), effective sterilization may be achieved by using the following parameters:

### European Sterilization Parameters

For users in Europe, Orthofix recommends the following parameters:

For Firebird Spinal Fixation System Cases 44-9010, 44-9011, 44-9012, 44-9013, 44-9020, 44-9030, 44-9040, 44-9050 and 61-9060, The NewBridge Laminoplasty Fixation System Instruments, NGage Surgical Mesh System Instruments, and ProView Minimal Access Portal (MAP) System Instruments, or when risk of prion contamination is suspected\*:

|                             |            |
|-----------------------------|------------|
| <b>Sterilization Method</b> | Steam      |
| <b>Cycle</b>                | Pre-vacuum |
| <b>Temperature</b>          | 134°C      |
| <b>Exposure Time</b>        | 18 minutes |
| <b>Minimum Drying Time</b>  | 30 minutes |
| <b>Number of Pulses</b>     | 4          |

\*Based on WHO (World Health Organization) guidelines

For all other instruments:

|                             |            |
|-----------------------------|------------|
| <b>Sterilization Method</b> | Steam      |
| <b>Cycle</b>                | Pre-vacuum |
| <b>Temperature</b>          | 134°C      |
| <b>Exposure Time</b>        | 3 minutes  |
| <b>Minimum Drying Time</b>  | 30 minutes |
| <b>Number of Pulses</b>     | 4          |

For all other implants:

|                             |            |
|-----------------------------|------------|
| <b>Sterilization Method</b> | Steam      |
| <b>Cycle</b>                | Pre-vacuum |
| <b>Temperature</b>          | 134°C      |
| <b>Exposure Time</b>        | 3 minutes  |
| <b>Minimum Drying Time</b>  | 30 minutes |
| <b>Number of Pulses</b>     | 4          |

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

### Non-European Sterilization Parameters:

For users outside of Europe, use one of the following recommended cycles with a FDA cleared sterilization wrap.

For Firebird Spinal Fixation System Cases 44-9010, 44-9011, 44-9012, 44-9013, 44-9020, 44-9030, 44-9040, 44-9050 and 61-9060, The NewBridge Laminoplasty Fixation System Instruments, NGage Surgical Mesh System Instruments, and ProView Minimal Access Portal (MAP) System Instruments:

|                             |               |
|-----------------------------|---------------|
| <b>Sterilization Method</b> | Steam         |
| <b>Cycle</b>                | Gravity       |
| <b>Temperature</b>          | 250°F (121°C) |
| <b>Exposure time</b>        | 30 minutes    |

or:

|                             |               |
|-----------------------------|---------------|
| <b>Sterilization Method</b> | Steam         |
| <b>Cycle</b>                | Prevac        |
| <b>Temperature</b>          | 270°F (132°C) |
| <b>Exposure time</b>        | 8 minutes     |

For all other instruments:

|                             |               |
|-----------------------------|---------------|
| <b>Sterilization Method</b> | Steam         |
| <b>Cycle</b>                | Gravity       |
| <b>Temperature</b>          | 270°F (132°C) |
| <b>Exposure time</b>        | 15 minutes    |
| <b>Drying time</b>          | 30 minutes    |
| <b>Double wrapped</b>       |               |

or:

|                             |                             |
|-----------------------------|-----------------------------|
| <b>Sterilization Method</b> | Steam                       |
| <b>Cycle</b>                | Prevac                      |
| <b>Preconditioning</b>      | Per manufacturer's settings |
| <b>Temperature</b>          | 270°F (132°C)               |
| <b>Exposure time</b>        | 4 minutes                   |
| <b>Drying time</b>          | 30 minutes                  |
| <b>Double wrapped</b>       |                             |

For all other Implants:

|                             |               |
|-----------------------------|---------------|
| <b>Sterilization Method</b> | Steam         |
| <b>Cycle</b>                | Gravity       |
| <b>Temperature</b>          | 270°F (132°C) |
| <b>Exposure time</b>        | 15 minutes    |
| <b>Drying time</b>          | 30 minutes    |
| <b>Double wrapped</b>       |               |

or:

|                             |                             |
|-----------------------------|-----------------------------|
| <b>Sterilization Method</b> | Steam                       |
| <b>Cycle</b>                | Prevac                      |
| <b>Preconditioning</b>      | Per manufacturer's settings |
| <b>Temperature</b>          | 270°F (132°C)               |
| <b>Exposure time</b>        | 4 minutes                   |
| <b>Drying time</b>          | 30 minutes                  |
| <b>Double wrapped</b>       |                             |

Sterilization in Rigid Sterilization Containers:

When using rigid sterilization containers, clean, inspect and prepare the rigid sterilization container according to the manufacturer's instructions. Select the appropriate rigid sterilization container with either filtered bottom or solid bottom to properly enclose the Orthofix case(s) (recommended 23¼" long x 11¼" wide container).

The following sterilization cycle has been validated:

|                             |                             |
|-----------------------------|-----------------------------|
| <b>Sterilization Method</b> | Steam                       |
| <b>Cycle</b>                | Prevac                      |
| <b>Temperature</b>          | 270°F (132°C)               |
| <b>Preconditioning</b>      | Per manufacturer's settings |
| <b>Exposure time</b>        | 4 minutes                   |
| <b>Drying time</b>          | 30 minutes                  |

## Non-Sterile Surgical Instruments and Implants Reprocessing Guide

### Packaging Integrity:

If the case, tray, rigid container, or sterile wrap is damaged, the product should not be used.

The product must be handled, stored, and opened in such a way that it is protected from inadvertent damage or contamination. If a loaner or consignment system is used, all sets should be carefully checked for completeness, and all components should be carefully checked for damage before use.

### Operative Technique:

Recommended operative techniques are available upon request from Orthofix at the phone numbers provided in the manufacturer's contact information above.

### Product Complaints:

Any Health Care Professional (e.g., customer or user of this system of products) who has any complaints or who has experienced any dissatisfaction with the product quality, identity, durability, reliability, safety, effectiveness and/or performance, should notify Orthofix Inc., 3451 Plano Parkway, Lewisville, TX 75056, USA, by telephone at 1-214-937-3199 or 1-888-298-5700 or by e-mail at [complaints@orthofix.com](mailto:complaints@orthofix.com).

### Notice about Serious Incidents:

Report any serious incident involving a device to Orthofix and the appropriate governing body in which the user and/or patient is established.

**NOTE:** Send any unwanted device back to Orthofix in a leak-proof container with the date, surgeon, and any reason for return.

### Latex Information:

The implants, instruments and/or packaging material for the [Device System Name] are not formulated with and do not contain natural rubber. The term "natural rubber" includes natural rubber latex, dry natural rubber, and synthetic latex or synthetic rubber that contains natural rubber in its formulation.

|                |                                                                                      |                                                                                                  |
|----------------|--------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|
| <b>Rx Only</b> | Federal (U.S.A.) law restricts this device to sale by or on the order of a physician |                                                                                                  |
| <b>EU REP</b>  | Authorized Representative                                                            |  Manufacturer |

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