

# NORTHSTAR™ FACET FUSION

FACET PREP INSTRUMENTS & OSTEOSTRAND® CERVICAL PELLETS  
SALES BROCHURE

**NorthStar™ Facet Fusion** is a novel procedural solution for reproducible biologic placement in the facet joint and lateral mass screw pathway preparation. Facet Prep Instruments and OsteoStrand® Cervical Pellets are packaged together within one product box to deliver a complete fusion kit, intended for use in conjunction with posterior cervical fixation systems. Featuring a set of single-use instruments and precision-shaped demineralized bone fiber pellets, NorthStar Facet Fusion enables:



**ACCURATE & CONTROLLED**

lateral mass screw preparation



**MAXIMIZED**

decorticated surface area  
within the facet joint




**SIMPLIFIED**

biologic delivery with precision-  
shaped demineralized bone  
fiber pellets



**EFFICIENT & REPRODUCIBLE**


method for robust facet fusion

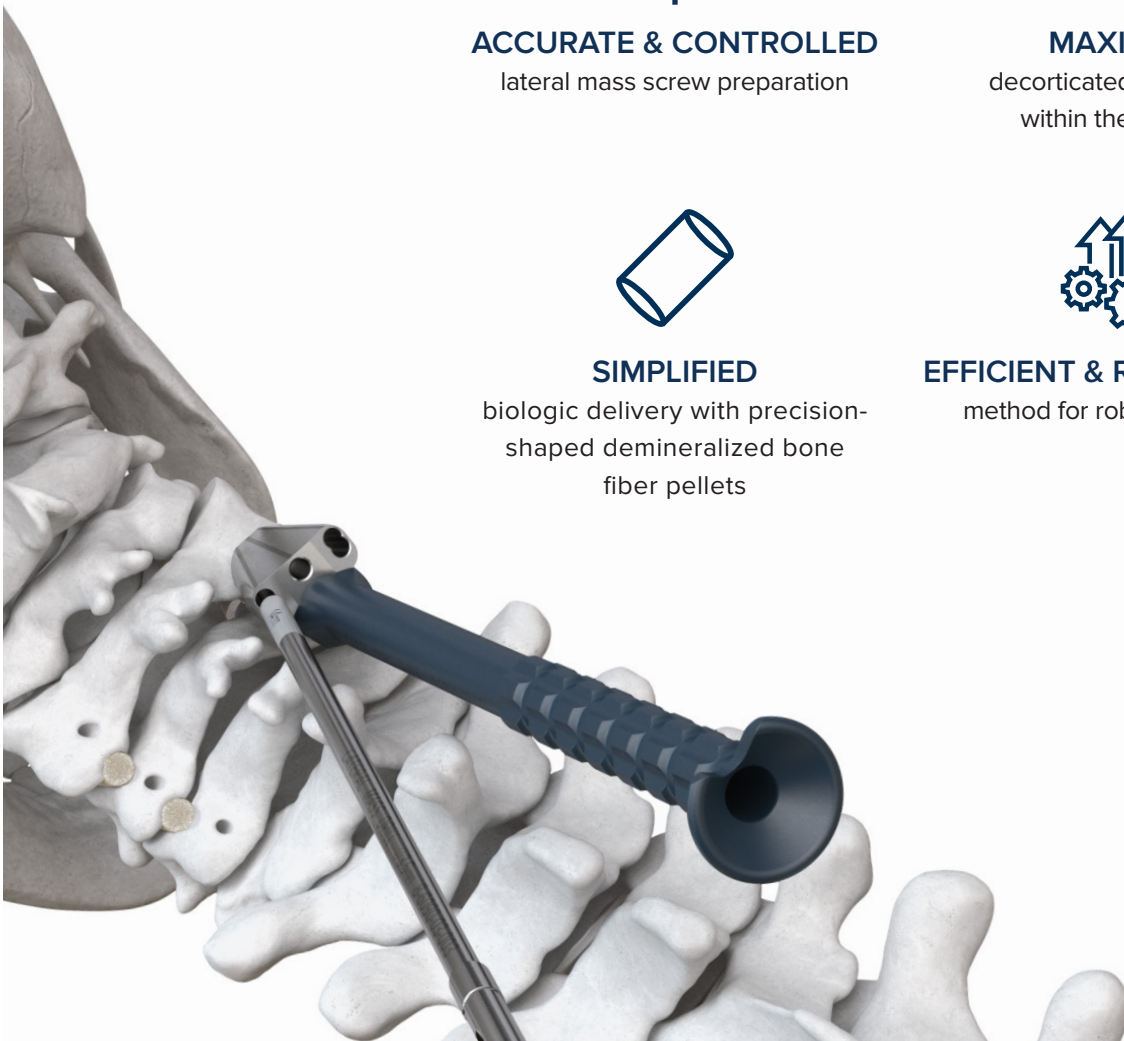

# NORTHSTAR™ FACET FUSION

FACET PREP INSTRUMENTS & OSTEOSTRAND® CERVICAL PELLETS

## FACET PREP INSTRUMENTS

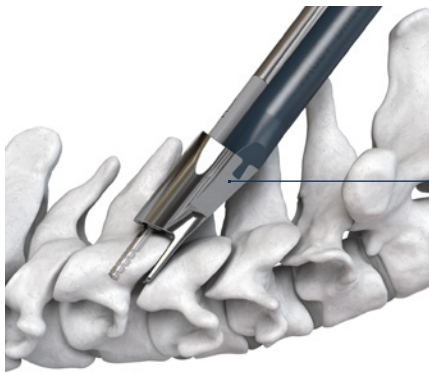

### Cannulated Facet Drill Guide

accesses the facet joint and positions cervical screw pilot holes




### 5.0mm Facet Drill

creates a precise defect for facet joint preparation

### Lateral Mass Drill

(2.25mm) for 3.5mm screws, with adjustable Depth Sleeves

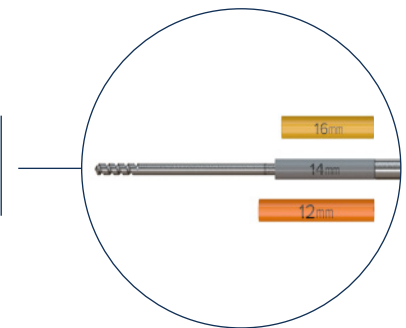


### Bone Graft Tamp

delivers biologic into the facet joint

## OSTEOSTRAND® CERVICAL PELLETS

- 100% demineralized bone fibers for maximized osteoinductive potential
- Entangled fibers provide osteoconductive scaffold
- Gentle expansion upon hydration *in situ* enables complete facet defect fill for fusion
- Available in 4- or 10-packs

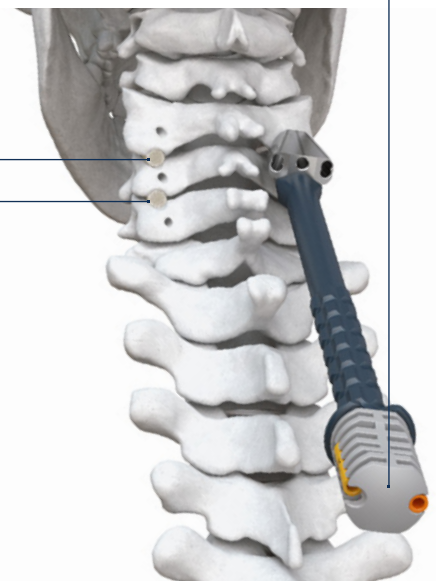


# NORTHSTAR™ FACET FUSION


FACET PREP INSTRUMENTS & OSTEOSTRAND® CERVICAL PELLETS

Compatibility with SeaSpine®  
posterior cervical  
fixation system

## NorthStar® OCT

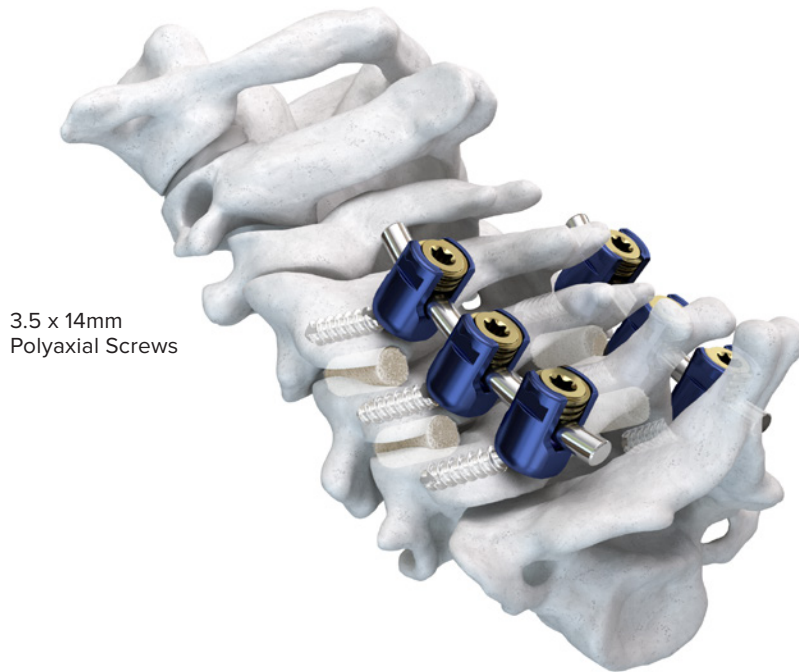

3.5 x 14mm  
Polyaxial Screws

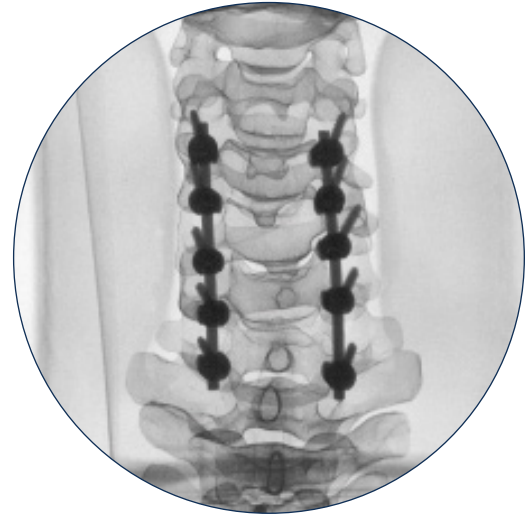

OsteoStrand® Cervical Pellets replaced with  
metal dowels for demonstration purposes



For more information or to place an order, please contact:  
TEL 866.942.8698 | FAX 877.558.6227  
customerservice@seaspine.com | [seaspine.com](http://seaspine.com)

#### Outside USA

TEL +1.760.727.8399 | FAX +1.760.727.8809  
INTERNATIONAL INQUIRIES [intlcustomer@seaspine.com](mailto:intlcustomer@seaspine.com)  
ALL OTHER INQUIRIES [customerservice@seaspine.com](mailto:customerservice@seaspine.com)

**IsoTis OrthoBiologics, Inc.**  
2 Goodyear, Irvine CA 92618  
TEL 800.550.7155 | FAX 800.471.3248 | WEB [seaspine.com](http://seaspine.com)  
IsoTis OrthoBiologics, Inc. is a member of the SeaSpine  
Orthopedics Corporation family of companies.  
Made in the U.S.A.



Warning: Applicable laws restrict these products to sale by or on the order of a physician. SeaSpine, the SeaSpine logo, IsoTis, NorthStar and OsteoStrand are registered trademarks of SeaSpine Orthopedics Corporation or its subsidiaries in the United States and/or other countries. NorthStar is a trademark of SeaSpine Orthopedics Corporation or its subsidiaries in the United States and/or other countries. ©2021 SeaSpine Orthopedics Corporation. All rights reserved. D0004137B 01