


# REEF<sup>®</sup> TA

TRANSFORAMINAL LUMBAR INTERBODY  
SALES BROCHURE

# REEF TOPOGRAPHY<sup>®</sup>



Undercut macrostructures designed to promote bony interlocking<sup>5,+</sup>

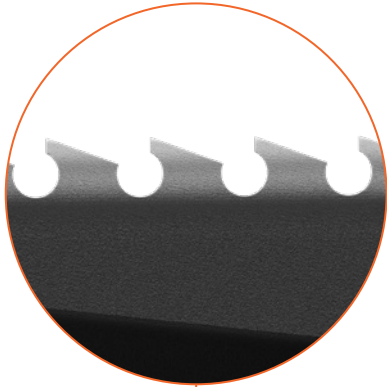

## UP TO 28% MORE

### NanoMetalene<sup>®</sup> Endplate Surface Area\*

Endplate features place graft material in direct contact with the endplates.

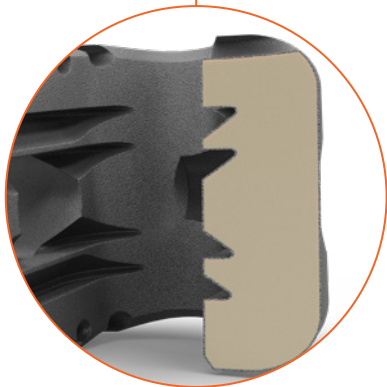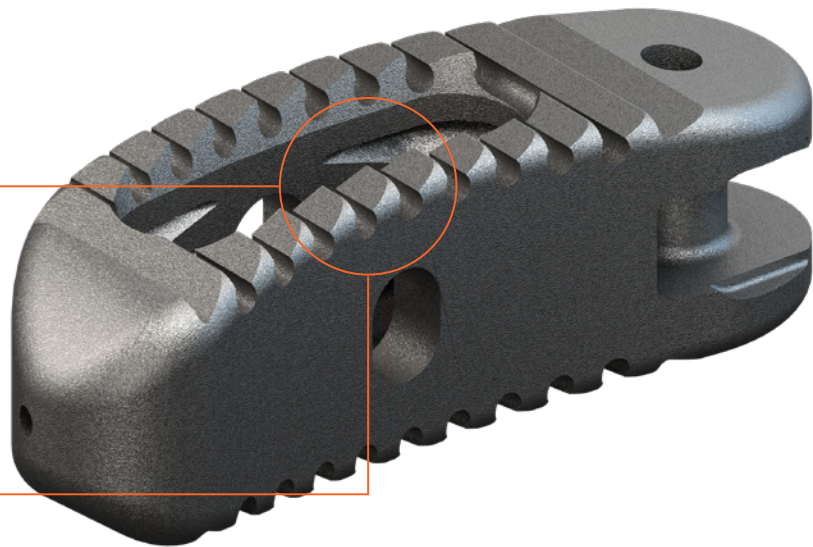

## UP TO 66% MORE

### NanoMetalene Aperture Surface Area\*


Aperture features secure graft within the aperture during interbody placement.

\*Data on file.


# NANOMETALENE<sup>®</sup> TECHNOLOGY

## Submicron titanium layer molecularly bonded to entire PEEK implant

Titanium surfacing resists wear debris.<sup>1</sup>


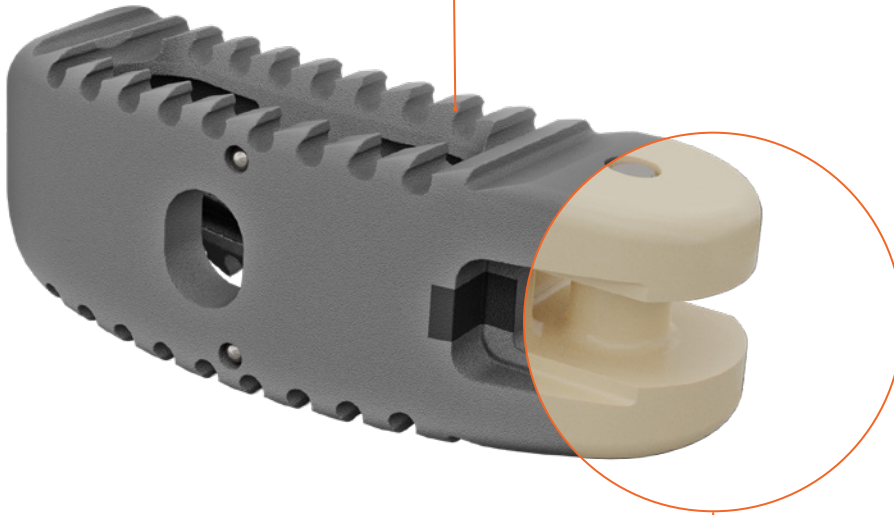


Mechanical properties of PEEK unaltered, providing stiffness on par with bone.<sup>1</sup>

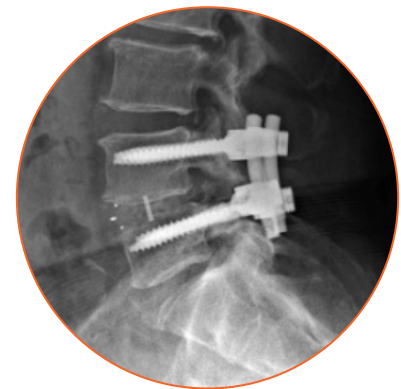

Radiolucent for post-op fusion assessment.<sup>2</sup>

## Preclinical results show greater bone ongrowth on NanoMetalene<sup>®</sup> vs. PEEK<sup>3,†</sup>

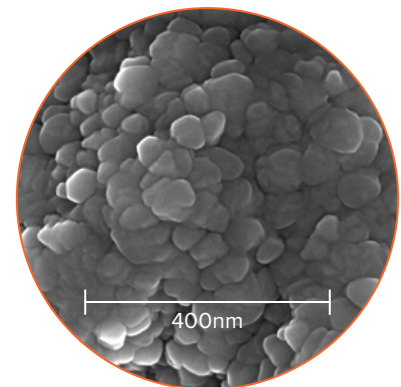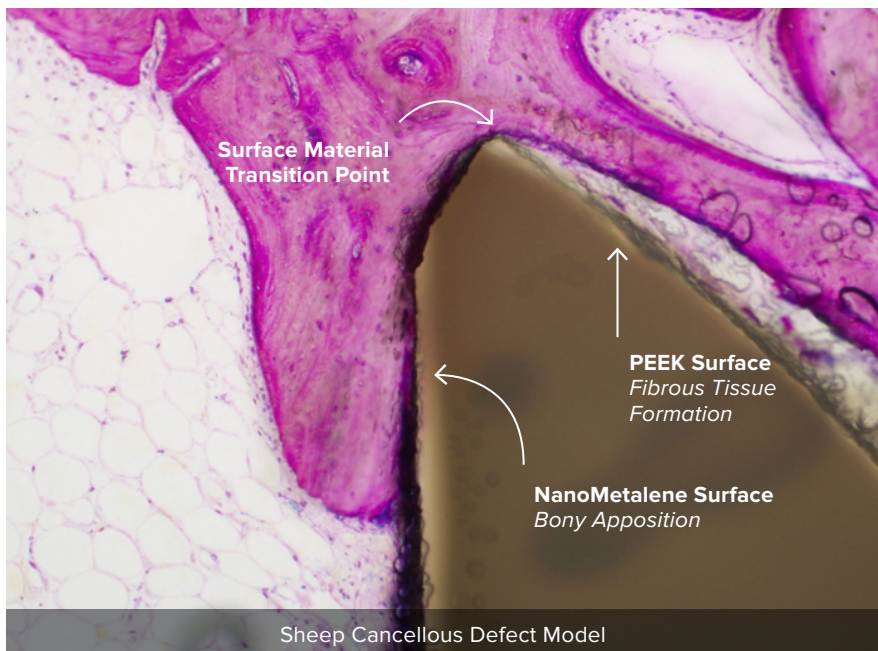

Rough submicron topography encourages integration.<sup>3,4</sup>

# REEF® TA

The SeaSpine® Reef® TA (TLIF Articulating) System is designed to reliably deliver an interbody to the anterior cortex during transforaminal lumbar interbody procedures. The implant and instrumentation provides a robust connection for implant delivery and articulation within the disc space until desired placement is reached. Implants feature NanoMetalene® surface technology and Reef Topography® for greater titanium surface area and improved biomechanical stability while maintaining the physical properties of PEEK.<sup>1,2,†</sup>

## Interbody Features

- Reef Topography and NanoMetalene surface technology
- Footprint options: 9 x 28mm, 9 x 32mm, 11 x 36mm
- Height options: 8–16mm\*
- Lordotic options: 8°, 15°



## Instrumentation to allow for up to 70° articulation



\*Note: Min./max. heights vary by footprint.




# PRECLINICAL EVALUATION



## Preclinical Evaluation of Reef Topography®

Clinically relevant endplate-sparing sheep interbody fusion model results comparing NanoMetalene® (NM) implants with and without Reef Topography.

### Endplate Undercut Macrostructures

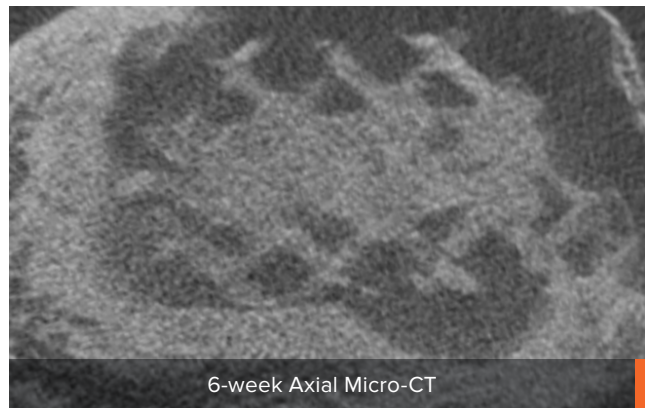


### Aperture Undercut Macrostructures

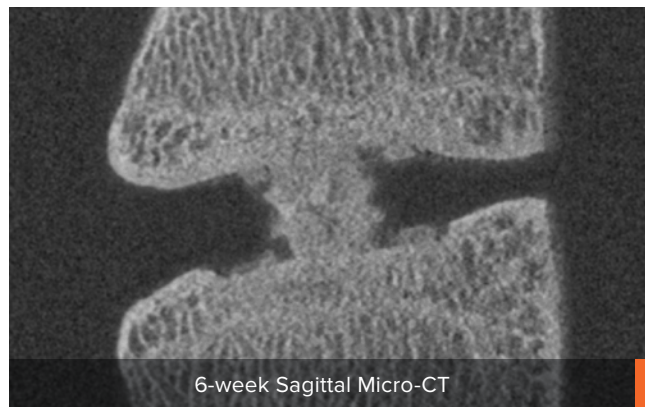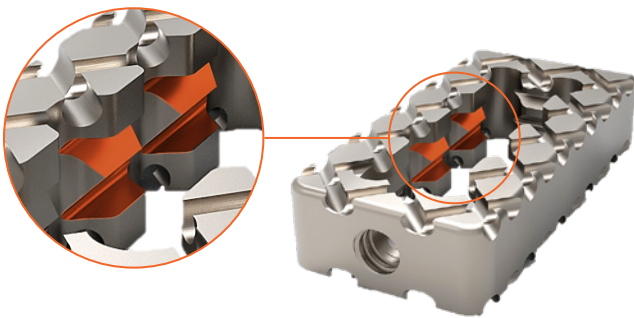

Undercut Macrostructures results in

**~3x increase in mechanical stability<sup>5,†</sup>**

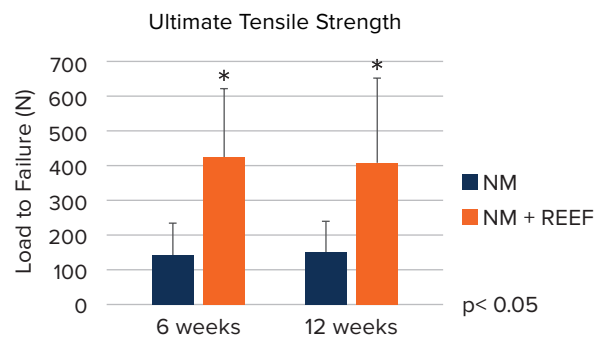

REEF TOPOGRAPHY MAY DRIVE EARLIER AND IMPROVED BIOMECHANICAL STABILITY<sup>5,†</sup>

# FUSION ENGINEERED®

## Deliberate Design. Driven by Science.

SeaSpine® is dedicated to providing a comprehensive portfolio of innovative, procedurally-focused products strategically designed to work together to drive fusion. The latest advancements in bone biology and materials science guide the development of in-house manufactured advanced orthobiologics and proprietary spinal implant technology engineered to address the many nuances of spinal pathology. SeaSpine products can be tailored to meet individual patient needs, delivering both clinical and economic value to patients, surgeons, and hospital systems.



### COMPLEMENTARY TECHNOLOGIES

The combination of SeaSpine spinal implant technologies and orthobiologics capitalizes on both design and science, unleashing a favorable environment for fusion.



### PROCEDURAL VERSATILITY

The comprehensive and procedurally-focused SeaSpine product portfolio offers flexibility that allows a tailored approach to each unique surgical situation.



### SCIENTIFICALLY DRIVEN

SeaSpine is committed to relying on data to drive each facet of our technologies toward the ultimate goal of fusion.

\*Preclinical testing, such as animal studies, may not be indicative of human results.

†Results from mechanical testing. Data on file. TR-0010-11-01


‡Results from imaging study. Data on file. TR-0010-11-01

§Walsh, et al. The in vivo response to a novel Ti coating compared with polyether ether ketone: evaluation of the periphery and inner surfaces of an implant. Spine Journal 2018 Jul; 18(7): 1231-1240

¶NanoMetalene SEM images on file. TR-0094-19-01

‡Results from preclinical in vivo testing. Data on File. D0003269

For more information or to place an order, please contact:  
TEL 866.942.8698 | FAX 877.558.6227  
customerservice@seaspine.com | [seaspine.com](http://seaspine.com)

 **SeaSpine Orthopedics Corporation**  
5770 Armada Drive  
Carlsbad, CA 92008 USA  
TEL 760.727.8399 USA | FAX 760.727.8809



#### Outside USA

TEL +1.760.727.8399 | FAX +1.760.727.8809

INTERNATIONAL INQUIRIES [intlcustomer@seaspine.com](mailto:intlcustomer@seaspine.com)

ALL OTHER INQUIRIES [customerservice@seaspine.com](mailto:customerservice@seaspine.com)